

High Resolution 4K Camera

MKC-750UHD

4%





4K RESOLUTION

4K High Resolution

(Unit:pixel)

By equipping the CCU with 4K (3840x2160) digital processing, the MKC-750UHD produces superb resolution with four times the sample rate of conventional HD (1920x1080).

(Horizontal Resolution 1600TVL, and S/N 58dB)

Quarter (1/4) video capture for full HD (1920x1080) display





ABOUT 4K

 $4 \rm K~(3840x2160)$ delivers four times as much detail as Full HD (1920x1080), providing high quality images as real as seen by the human eye. Especially on a big screen, it reproduces images superior to Full HD.





COMPACT AND LIGHTWEIGHT

Compact and Lightweight Camera Hea

Although the MKC-750UHD is a 4K camera, it has a very compact and lightweight camera head.

FEATURES & FUNCTIONS

High Sensitivity

When providing 4K resolution, the MKC-750UHD has a sensitivity of F11 (LINE MIX ON) at 2000 lx.

(Various Image Outputs)

The CCU outputs digital HD signals (HD-SDI and DVI), as well as, a 4K signal.

Automatic Iris Control

By connecting an exclusive camera adaptor, the auto iris control function can be used. Optical Iris in this connection.

Automatic Electronic Shutter System (Shutter speed)

The continuously variable electronic shutter can be used with a manual iris lens. This Iris functionality maintains optimized performance over a wide range of lighting conditions.

Auto Gain Control

The gain can be adjusted automatically to provide the proper picture brightness when the light level is not enough.

Flip and Turn Function

The images can be flipped horizontally and vertically by the digital processing to correct the TV image against the inversion caused by optical instruments. Remote control is available via an optional external unit such as a foot switch.





Antibacterial Coating

The camera head and CCU (camera control unit) are coated to improve infection control.

Picture Freeze

By pressing the Freeze button on the CCU, it is possible to hold the current image.

Scene File

4 scene files are provided to keep pre-set memories of various control functions for convenient use. Scene file data can be saved into and loaded from a USB memory.

I Function Switch

Four switches are available on the CCU to custom assign to control functions important to the user.



Digital Zoom

Electronic zoom for the center of an image is available, up to four times.

Level Indicator

Video level can be checked and adjusted using a level indicator on the CCU.



ACCESSORIES

Surgery Microscope Adaptor

CFA-400R
Focus/Remote control



CFA-400



Other Accessories

Camera cable

Single: MCCGHPP-050 (5m) Single: MCCGHPP-070 (7m)

Combination: MCCGHPP-050 (5m) +

MCCGCJP-100 (for extension 10m)





Specifications

Rating		
Optics	RGB prism	
Image Pickup Device	1/3-inch, CMOS x 3	
Lens Mount	C-Mount	
Output Signals	4K	3G-SDI (Level A) × 4 BNC 1ch
		3840 × 2160p
	HDTV	HD-SDI 2ch
		1920 × 1080i/p
		DVI 1ch
		1920 × 1080i/p
External Reference	HDTV Tri-Sync 1ch	
Camera cable (option)	Single: MCCGHPP-050 (5m)	
	Single: MCCGHPP-070 (7m)	
	Combination : MCCGHPP-050 (5m) + MCCGCJP-100 (for extension 10m)	
Dimensions / Weight	Camera Head	W34 × H40 × D40mm (excluding protrusions)
		/ 100g approx.
	Camera Control Unit	W210 × H80 × D270mm (excluding protrusions)
		/ 3kg or less
Power Requrerment	AC 100V - 230V ±10% 50/60 Hz	
Power Consumption	29W	
Safety Approval	EN60601-1, ANSI/AAMI ES60601-1	

Performance	
Resolution	1600 TV Line
S/N	58dB (Gamma/Detail OFF)
Sensitivity	2000 lx F11 /3200K or more (LINE MIX ON)

Function		
Electronic Shutter	Manual 1/100~1/10000	
(Rolling Shutter)	Auto 1/100~1/10000	
	Slow 1/4, 1/8, 1/15, 1/30 field/sec	
Gain	-9dB~24dB (3dB steps)	
Automatic Functions	Auto White Ballance Auto Gain Control Auto Iris Cointrol Auto knee	
Image Correction Functions	Shading Correction, DTL Correction, Flare Correction, Auto knee	
Auto Iris Sample Area	Narrow / Middle / Full / Circle / USER	
Image Flip and Turn	Horizontal, Vertical	
Scene File	4 files (by CCU front switch) Setting data can be stored in USB memory	
H/V Centering	Built-in	
Digital Zoom	1.1x -4x	
Color Bar	Built-in	
Color Correction (paint)	R-channnel / B-channel Correction available via CCU front switch	
Brightness Adjustment	Available via CCU front switch	
Freeze Function	Built-in (Controlled by CCU front switch or optional foot switch),	
	still image / live image can be exchanged	
Image Level Indicator	7 steps (by LED display on CCU)	

Front & Rear Panel


Front Panel




Rear Panel





Dimensions (Unit:mm)

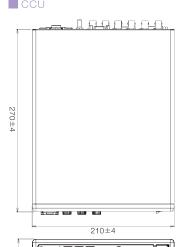


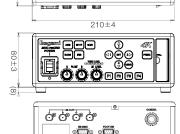

૽ ૽૽૽૽૽૽૽ઌ૿૽ઌ૿૽ઌ૿૽૽ૼ

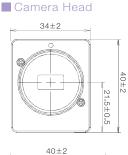

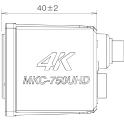

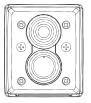

Design and specifications are subject to change without notice.



Ikegami Tsushinki Co.,Ltd.

■ URL: http://www.ikegami.co.jp/en/

U362A168-EU

■ Head Office : 5-6-16 Ikegami, Ohta-ku, Tokyo 146-8567, Japan · · · · · TEL.+81-(0)3-5700-1111 / FAX.+81-(0)3-5700-1117 ■ Overseas Sales Division: 5-6-16 Ikegami, Ohta-ku, Tokyo 146-8567, Japan · · · · · TEL.+81-(0)3-5700-4117 / FAX.+81-(0)3-5748-2200

■ Beijing Office: Room No.307 Tower B, Jia Tai International Mansion No.41East Fourth Ring Road, Chaoyang District, Beijing China ···TEL.+86-(0)10-85710350 / FAX.+86-(0)10-85711804





Ikegami Electronics (U.S.A.), INC.

Ikegami Electronics (Europe) GmbH

■ URL: http://www.ikegami.de

■ URL: http://www.ikegami.com

......TEL.+49-(0)2131-1230 / FAX.+49-(0)2131-102820



■ Headquarters : Ikegamistraße 1, D-41460 Neuss, Germany

Ikegami Electronics Asia Pacific Pte.Ltd. ■ URL: http://sg.ikegami.co.jp ■ Headquarters: 1 Tampines Central 5, #06-04 CPF Tampines Building, Singapore 529508 · · · · · · · · · · · · TEL.+65-6260-8820 / FAX.+65-6260-8896